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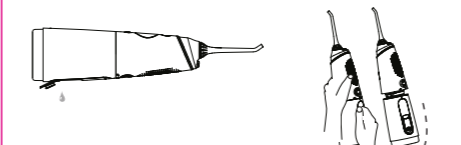
**11. TROUBLESHOOTING**

PROBLEM	PROBABLE CAUSE	SOLUTION
Water leaking between the jet tip end and the irrigator handle	The jet tip is not inserted correctly	Pull the tip out and connect it to the socket
Water tube is leaking	Damaged water tube	Contact the seller or the manufacturer service center
Water leaking from the tank	No valve at the end of the tank, or the tank is installed the wrong way around	Mount the valve in the correct way or report the problems to the seller
Water pressure is very low	The tip is clogged or deformed	Replace the tip
Irrigator is not turning on	The device is damaged, or the battery is flat	Check the battery charge level or contact the seller
Fracture of the gum jet tip	The jet tip can be used	Replace the jet tip for a new one
The buttons on the irrigator handle are not working	The battery is flat or installed incorrectly	Check or change the batteries

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
**Importer:** KGK Trade sp. z o.o. sp. k. ul. Ujastek 5B 31-752 Kraków Polska

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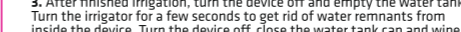**Dear Customer,** Thank you for trusting our brand and buying Evorei Boost 1, an irrigator used for oral hygiene and care. The device is not a medical product and is designed for household use only.

Before using the device, please follow the user manual carefully. Failure to comply with the rules presented in the instructions can cause injuries or device damage. The manufacturer is not liable for the damages resulting from not following the instructions.

**01. SAFETY**

1. Children of 8 years and older and people with reduced motor, mental or cognitive abilities can use the device only under guardian supervision.
2. It is not recommended to use the device longer than 5 minutes or more often than every 2 hours.
3. Do not use the device if you have sores or cuts inside your mouth.
4. Remove all jewelry from your mouth before using the device.
5. Do not toss or put any objects into the device when open.
6. If your doctor recommended you take antibiotics before doing any dental work, consult your dentist before using the device.
7. The device is not a toy.
8. Cleaning and maintenance should not be performed by children without adult supervision.
9. Do not disassemble the device except for when it is supposed to be utilized.
10. Do not insert the plug into the socket with a different voltage than the one specified for this device.
11. Do not fill the tank with water over 40°C!
12. Check the device before every use. If it appears to be damaged, do not use it.
13. Do not hit or drop the device. If you drop or hit the device, and it becomes damaged – stop using it immediately.
14. Store the charging cord away from heated surfaces.
15. Do not power the device without previously mounting a jet tip.
16. Do not immerse the device in water or other liquids.
17. Keep the tip of the water flosser out of reach of children – choking hazard!
18. Mouthwash can damage the device. Do not pour mouthwash inside the water flosser tank.
19. If you feel pain or discomfort during or after irrigation, or experience any other reactions to the water flosser, consult your doctor.
20. Do not press the tip of the flosser too hard against your gums or teeth.
21. Use only original tips included in the set.

**22.** The device is designed to clean the spaces between teeth and tongue and to improve gum health. Do not use the device for other reasons than what stated in the user manual!

**23.** Do not direct the water stream directly under the tongue, inside the ear, nose, and other delicate places. The device can produce such pressure that can damage these sensitive spots.

**24.** The device is not designed to be used by children and people with reduced mental or cognitive capabilities.

**25.** Do not try and fix the device on your own or modify and disassemble the device parts.

**26.** Do not put the device in a place where it can fall into a tub or sink.

**02. ORAL HYGIENE IRRIGATOR**

Designed to clean the spaces between permanent braces, crowns, and dental bridges. It gets rid of up to 99% plaque and perfectly cleans the spaces between teeth and below gum lines.

**Standard teeth cleaning:** Washing your teeth with a traditional brush does not provide a thorough removal of all the leftover food particles and bacteria from between the teeth.

The range of dental floss is also limited - it can reach only the upper gum line.


Irrigator will precisely remove food particles, residue, and pathogenic bacteria from the spaces between teeth.

**03. INCLUDED**


1. 4x standard jet tips
2. 1x charging cable
3. 1x user manual

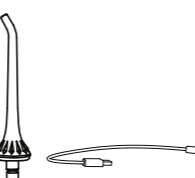**ATTENTION!** Some people may experience bleeding from gums when first using the Irrigator. Bleeding gums are a normal phenomenon which depends on the gum condition. If the bleeding does not stop after 2 weeks or regular irrigator usage, it can indicate that you are suffering from periodontal disease. Stop using the device and consult your doctor.

**04. DEVICE DESCRIPTION**


1. ON/OFF switch
2. 0-LED display
3. Mode button
4. Jet Tip
5. Lock button
6. Recharge input
7. Water tank

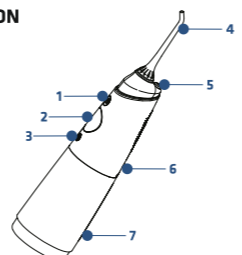**05. BEFORE USE**

**IMPORTANT!** Irrigator will be fully ready to use after previous battery charge! Time of battery charge is around 4 hours. **DO NOT USE THE DEVICE WHEN CHARGING!**

1. Plug the charging USB cable into the power source.
2. When charging the battery, the display will show battery icon.
3. Device charging time is around 4 hours.

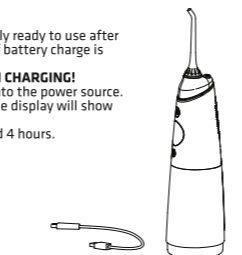**06. USING THE DEVICE**

1. Turn the device by pressing the ON/OFF button situated on the upper part of the case right above the screen.
2. After turning the device on, the irrigator will automatically switch to the last used settings.

**3.** The device is equipped with 5 irrigation modes (STRONG, MEDIUM, SOFT, PULSE and DIY). To change the working mode, press the button below the screen.

**4.** In DIY mode you can regulate the pressure level by pressing and holding the mode change button.

**5.** The device automatically turns the water stream off after working for 120 seconds in any mode.

**6.** Battery level can be checked on the device display. If the battery runs out of power, the display will not light up.

**7.** Before the first use and when device has not been used for a long time (at least 3 months), charge the device for 4 hours. Fully charged irrigator can work for around two weeks.

**07. IRRIGATION TIPS**

**1.** Standard irrigation tip: used to clean the spaces between teeth, gum gaps and gum pockets. Classic tip is designed for daily cleaning in hard-to-reach places.

**08. TEETH CLEANING TECHNIQUES**

1. Place the irrigator in a vertical position. Adjust the jet tip so the water stream is turned to a 90° angle to the gums to avoid splashing with water, and so the water can easily flow from your mouth to the sink.
2. To get the best results, start from the back section of the teeth and slowly move the jet tip towards the front. Move the tip along the gum line and stop shortly between the teeth lines. Continue the process until you fully clean the inner and outer sides of top and bottom teeth.