EVOREI

Sonic 2

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PL

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USERS MANUAL

EN


**Dear Customer,** thank you for purchasing Evorei Sonic 2. The device you purchased is a sonic toothbrush. **Before using the device, please read the user manual carefully.**

#### 01. SAFETY MEASURES

- If you have never used a sonic toothbrush before, you may experience slight discomfort and bleeding gums the first time you use the device. If these symptoms continue after two weeks of use, contact your dentist.
- Follow your dentist's recommendations concerning dental care.
- The device is used to brush teeth, gums and tongue. Do not use it for any other purpose. If you feel pain or discomfort while using the brush, stop using it and consult your dentist.
- Do not use the brush head if it is bent, distorted bristles. It is recommended to replace the brush head every 3 months. If the brush head is visibly damaged, replace it immediately with a new one.
- Do not throw the device or drop it. If the device is damaged, stop using it.
- Do not use the sonic toothbrush in the bath or shower.
- The product contains a built-in rechargeable battery. Do not throw the device into a fire and keep it away from high temperature areas. Do not charge or start the device in areas with high temperature.
- Use only the original accessories provided by the manufacturer.

#### 02. DEVICE DESCRIPTION

- Brush head
- Power/ mode change button
- Brush holder
- Mode indicator
- Battery status indicator
- Charging port (Behind the lid)
- Carrying case
- Charging cable

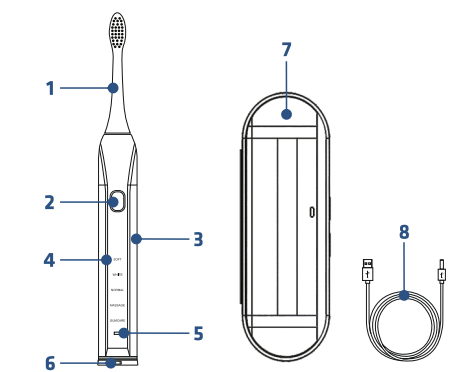
#### 03. PACKAGE CONTENTS

- Evorei Sonic 2 sonic toothbrush, Carrying case, 2x brush head, charging cable, user manual.

#### 04. CHARGING

- Uncover the charging port behind the lid at the bottom of the device.
- Connect the device to power using the charging cable.
- The battery status indicator will turn red while charging. When the device is fully charged, the indicator turns green.

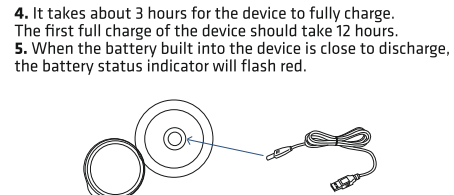
#### 05. FEATURES

- Power/ mode change button:** press to turn the device on or to change the mode of operation. Hold down the button for 5 seconds to turn off the device (the last selected mode will be selected when the device is turned on again). Pressing the button while GUMCARE mode is on will turn off the device.
- Operating modes:**
- SOFT** – gentle brushing, adapted for people who use the sonic toothbrush for the first time or who have delicate teeth.
  - WHITE** – whitening, the brush removes stains from tobacco, coffee and similar substances.
  - NORMAL** – standard brushing.
  - MASSAGE** – rhythmic massage of teeth and gums.
  - GUMCARE** – mode designed for people with missing teeth or infections of the mouth.
- Timer:** The device stops brushing every 30 seconds for 1 second to allow the user to move the brush to another part of their mouth. After 2 minutes, the device will turn off automatically.

#### 06. CLEANING AND MAINTENANCE OF THE DEVICE


- Cleaning the brush head:** wash under running water after each use and let the head dry naturally. Do not pull or press the bristles – this may distort them and make them fall out.
- Cleaning the device holder:** wipe the case with a delicate, wet cloth and let the device dry naturally. Do not submerge the device under water. Do not use water above 40°C or acidic cleaners.
- Replacement of the brush head:** it is recommended to replace the head every 3 months. Use only the heads recommended by the brush manufacturer.
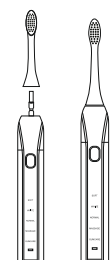
#### 07. USING THE DEVICE

- Place the brush head on the metal rod sticking out of the device holder. Make sure that the bristles are turned in the same direction as the power button.
- Squeeze the toothpaste onto the bristles and gently press them against your teeth at a 45° angle.

- Press the power button and start brushing your teeth. It is recommended to clean the teeth for 2 minutes. Every 30 seconds, the device stops vibrating for 1 second to allow the user to move the toothbrush to another part of the mouth.

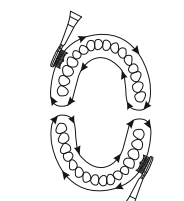- To avoid splashing toothpaste foam, turn off the toothbrush before pulling it out of your mouth.
- Rinse the brush head under running water and shake it several times to remove excess water. Wipe the handle of the device with delicate cloth.
- When you are not using the brush, keep it upright.

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#### 08. FAQ

- The brush does not work.** If the brush is new (used for less than 3 months), connect it to power and charge it for more than 24 hours. **Vibrations weaken during use.** You press the toothbrush too hard. It should gently touch your teeth. **After charging, the brush only works for a short while.** Try charging your device for more than 12 hours. If the problem persists, please contact the repair service.

#### 09. SPECIFICATION


- Vibration per program:** 78 000
- Operating modes:** Soft (gentle brushing), White (whitening brushing), Normal brushing, Massage, Gum care
- Number of modes:** 5
- Input voltage:** DC 5V - 1A
- Power:** 1 W
- Waterproof rating:** IPX7
- Battery:** 3.7 V 800 mAh, Li-ion
- Charging time:** Up to 3 h
- Working time:** 40-60 days (when washed twice a day)
- Brush weight:** 96 g
- Weight with case:** 190 g
- Brush dimensions:** 238x28x28 mm
- Case dimensions:** 205x65x25 mm

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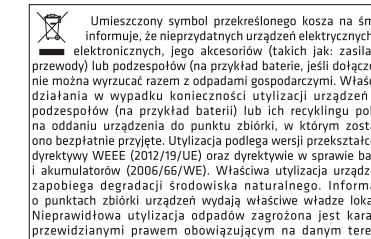**Producent:** KGK Trade sp. z o.o. sp. k. Os. Urocze 12 31-953 Kraków Polska

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
#### 10. WARRANTY/ CLAIMS

This product is covered by a 24-month warranty. The warranty conditions are available on the website: <https://evorei.pl/warranty/>

Complaints should be submitted with the complaint form which can be found at <http://reklamacje.kgktrade.pl/>

Contact details and the address of the service can be found on the website: [www.evorei.pl](http://www.evorei.pl)

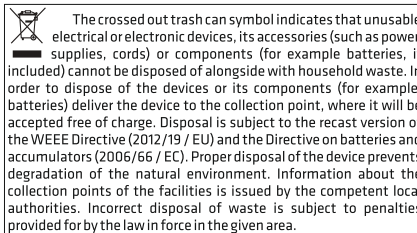**Manufacturer:** KGK Trade sp. z o.o. sp. k. Os. Urocze 12 31-953 Cracow Poland

Made in P.R.C.