[SONIC UV PRO]

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PL

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05. USING THE DEVICE

- Place the brush on the charger.
- When the brush is charging, the LED on the battery indicator on the brush will be red. When the device is fully charged, the LED turns green.
- When the battery of the toothbrush is close to being discharged, the battery indicator the device blinks with a red light.
- It takes 6 hours to fully charge the device. The first charge of the device should take 12 hours.

05. USING THE DEVICE

- Place the toothbrush head with the bristles on the metal rod sticking out of the handle. There should be a 1mm gap between the head and the body of the device.
- Squeeze toothpaste onto the toothbrush bristles.
- Press the brush gently to your teeth and gums.
- Press the power/mode switch to start cleaning your teeth. Clean every part of your teeth thoroughly.
- After brushing your teeth, wash the brush under running water. Shake the remaining water off the brush bristles and use a delicate cloth to wipe the handle of the device.
- When you are not using the brush, place it vertically upright on the wireless charger.

NOTE: to avoid splashing toothpaste foam around your bathroom, turn off the device before removing it from your mouth.

06. DEVICE MAINTENANCE

- CLEANING THE CASE**
 - Before cleaning the case, disconnect it from the power supply.
 - Using a delicate cloth, wipe and clean the case and charger. Pay special attention when cleaning the UV lamp.

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08. SPECIFICATION

Maximum quantity rotations per minute:	62000 SOFT – 62000 movements/min WHITE – 48000 movements/min NORMAL – 48000 movements/min MASSAGE – 42000 movements/min GUMCARE – 62000 movements/min
Charging time:	Up to 6 hours
Working time:	45-60 days (after washing twice during the day)
Battery:	3V, 800 mAh, lithium
Number of operating modes:	5
Water resistance:	IPX7
Brush weight:	100 g
Weight with case:	246 g
Brush dimensions:	238x28x28 mm
Dimensions with case:	218x77x31 mm

07. FAQ


- The brush does not work.**
 - If the brush is new (used for less than 3 months), connect it to power and charge it for more than 24 hours.
- Vibrations weaken during use.**
 - You press the toothbrush too hard. It should gently touch your teeth.
- After charging, the brush only works for a short while.**
 - Try charging your device for more than 12 hours. If the problem persists, please contact the repair service.

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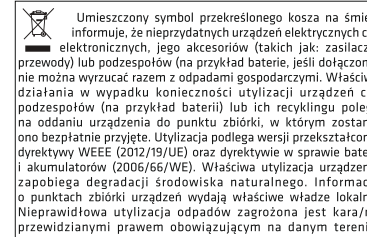USER'S MANUAL

EN

Dear Customer, thank you for purchasing **Xblitz EVORI SONIC UV PRO**. The product you have obtained is a sonic toothbrush with a UV-light disinfecting carrying case and a charging station. **Before using your device, please read the instructions below.**


01. SAFETY MEASURES

- If you do not understand the instructions for use or the precautions listed in this manual, do not use the device. If you do not use your device correctly, you may damage it.
- Do not disassemble the product or replace any of its components. This may result in product damage, electric shock and the warranty being voided.
- Do not use accessories other than those recommended by the manufacturer.
- Ultraviolet (UV) light can damage your eyes or skin. Do not expose any part of your body to the UV lamp placed in the device case.
- This product is for brushing teeth and gum massage only. Do not use the device for any other purpose.
- Do not connect the device to a voltage outlet other than one with parameters listed in the specification.
- Do not push the power plug into the power socket by force. If the plug does not easily enter the power socket, stop using the device.
- Stop using the device if it is clearly damaged.
- Keep the charging cord away from hot surfaces.
- If UV light is still on when the case is opened, stop using the case.
- Do not submerge the product under water. This can lead to a malfunction and may cause fire, electric shock or an explosion.

- Keep the device away from water containers such as a bathtub or sink.
- If the charger falls into the water, immediately disconnect the plug from the power outlet.
- Do not pull out the plug of the device when your hands are wet.
- Never use a charger if its cable or plug are damaged.
- Do not place the device in rooms with high moisture or high temperature, such as bathrooms.
- Disconnect the device from the power supply if you do not use it for an extended period of time to avoid the risk of fire.
- If the toothpaste used contains peroxides or bicarbonates (often present in whitening toothpastes), wash the device thoroughly after use to avoid destroying its plastic casing.
- If your gums bleed or feel uncomfortable for at least a week of you using the device, consult your dentist before further use.
- The device is not intended for the use by persons with limited physical, cognitive or mental capabilities or who lack experience or knowledge, including children. If such a person wishes to use the device, they must be supervised or instructed on the use of the device by the person responsible for them.
- The device is not a toy – keep it away from children.
- The device is used for personal hygiene and should not be used on many people in dentist's offices or other such professional environment.

02. DEVICE DESCRIPTION

- Power/ mode switch
- Mode indicator
- Battery indicator
- Toothbrush head
- Indicator of ongoing disinfection
- USB port
- Wireless charger (docking station)

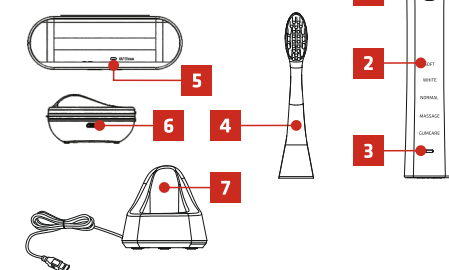03. PACKAGE CONTENTS

- Xblitz Evorei Sonic UV PRO
- Wireless charger (docking station)
- 4 toothbrush heads with bristles
- Carrying case for the device
- User's Manual

04. FUNCTIONS

1. BRUSH

- Switch on the device:** Press the power/ mode switch to turn on the device.
- Mode change:** Press the power/ mode switch to change the brush mode.

SOFT – soft, gentle brushing.
WHITE – teeth whitening.
NORMAL – normal brushing frequency.
MASSAGE – strong, massaging brushing.
GUMCARE – gums massaging.

After every 30 seconds of brushing, the brushing stops for 1 second to allow the user to change the position of the toothbrush. The device turns off automatically after 2 minutes, which is exactly for how long you should brush your teeth.

- Switch the device off:** Hold the power/ mode switch for 2 seconds in any mode to turn the device off.

2. UV CLEANING

- Connect the USB cable to the USB port on the case, connect the other end of the cable to USB power adapter (DC 5V).
- Close the case. The UV lamp should turn on automatically and the LED marked "UV Clean" should turn on, signaling that the disinfection of the brush is ongoing. The lamp will turn off automatically after 15 minutes, ending the disinfection. The UV disinfecting lamp is only turned on when the case is closed.

3. CHARGING


- Clean the toothbrush before charging.
- Connect the wireless charger to the power supply using a USB cable and a USB (DC 5V) adapter.

09. WARRANTY/ CLAIMS

This product is covered by a 24-month warranty. The warranty conditions are available on the website: <https://xblitz.pl/warranty/>

Complaints should be submitted with the complaint form which can be found at <http://reklamacje.kgktrade.pl/>

Contact details and the address of the service can be found on the website: www.xblitz.pl

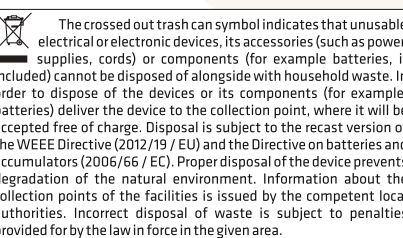Importer:
KGK Trade sp. z o.o. sp. k.
Os. Urocz 12
31-953 Cracow, Poland
Made in P.R.C.

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